


FLEXILINE PRO

FLEXILINE PRO

The most advanced design in a single use product!

- Flexible loops which may be shaped to the specific requirements of the patient
- Symmetrical loop for increased visibility of active tip
- Ultra-fine rigid tungsten wire for minimum eschar build up
- Non stick ball electrodes minimise adhesion.


Adjustable loops

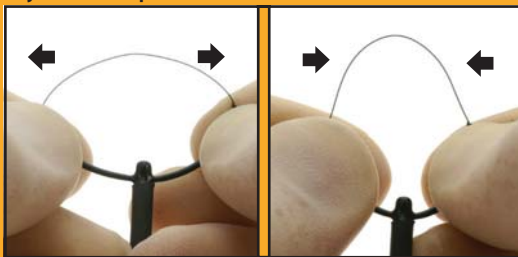

Non stick balls



Single use colposcopy loop electrodes for reduced adhesion and cleaner biopsies

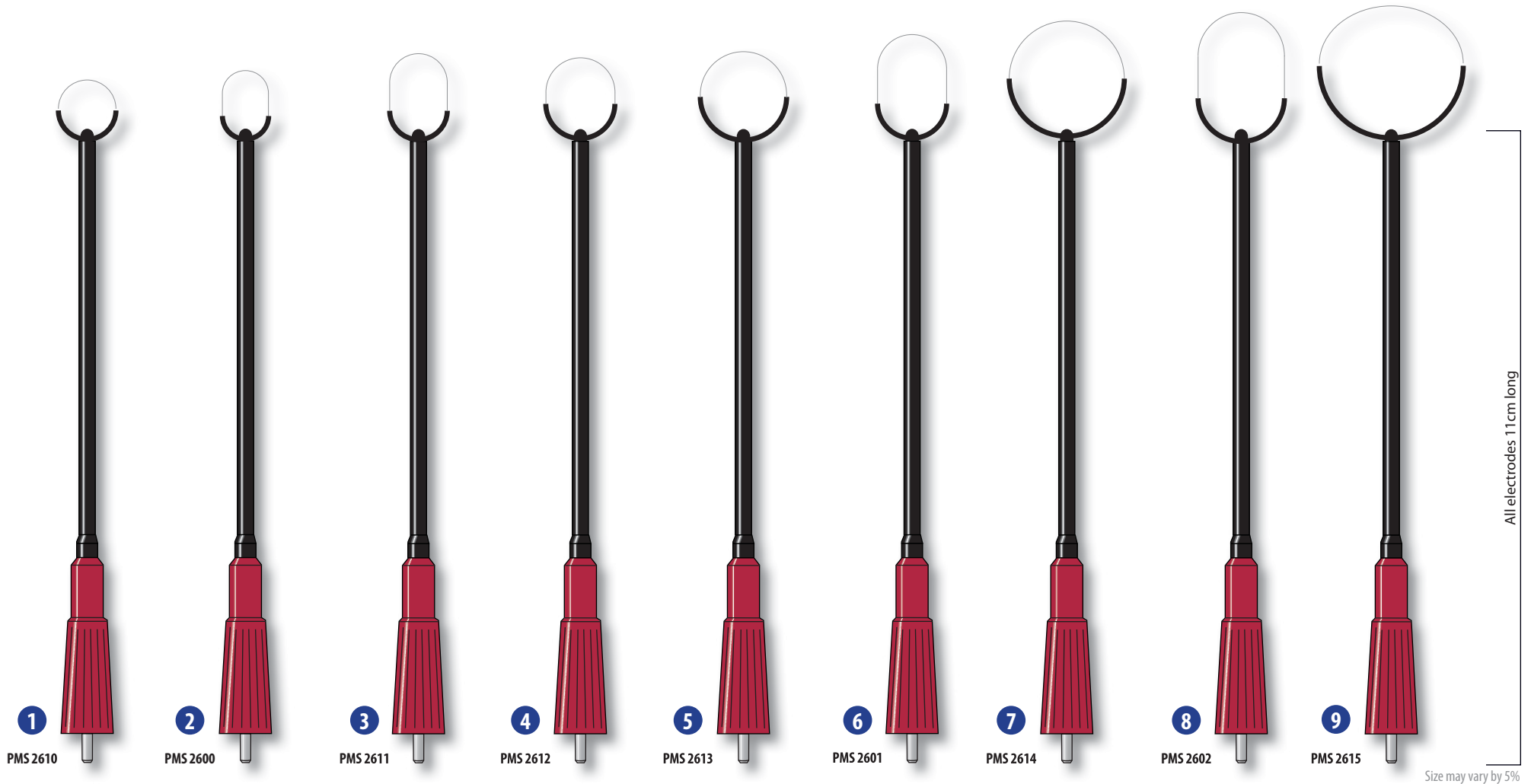


FLEXILINE PRO SINGLE USE LOOP ELECTRODES

- ① Loop electrode, 10 mm × 5 mm, 11 cm, pack 10 PMS2610
- ② Loop electrode, 8 mm × 8 mm, 11 cm, pack 10 PMS2600
- ③ Loop electrode, 10 mm × 10 mm, 11 cm, pack 10 PMS2611


- ④ Loop electrode, 12 mm × 8 mm, 11 cm, pack 10 PMS2612
- ⑤ Loop electrode, 15 mm × 8 mm, 11 cm, pack 10 PMS2613
- ⑥ Loop electrode, 12 mm × 12 mm, 11 cm, pack 10 PMS2601

- ⑦ Loop electrode, 20 mm × 10 mm, 11 cm, pack 10 PMS2614
- ⑧ Loop electrode, 15 mm × 15 mm, 11 cm, pack 10 PMS2602
- ⑨ Loop electrode, 25 mm × 10 mm, 11 cm, pack 10 PMS2615



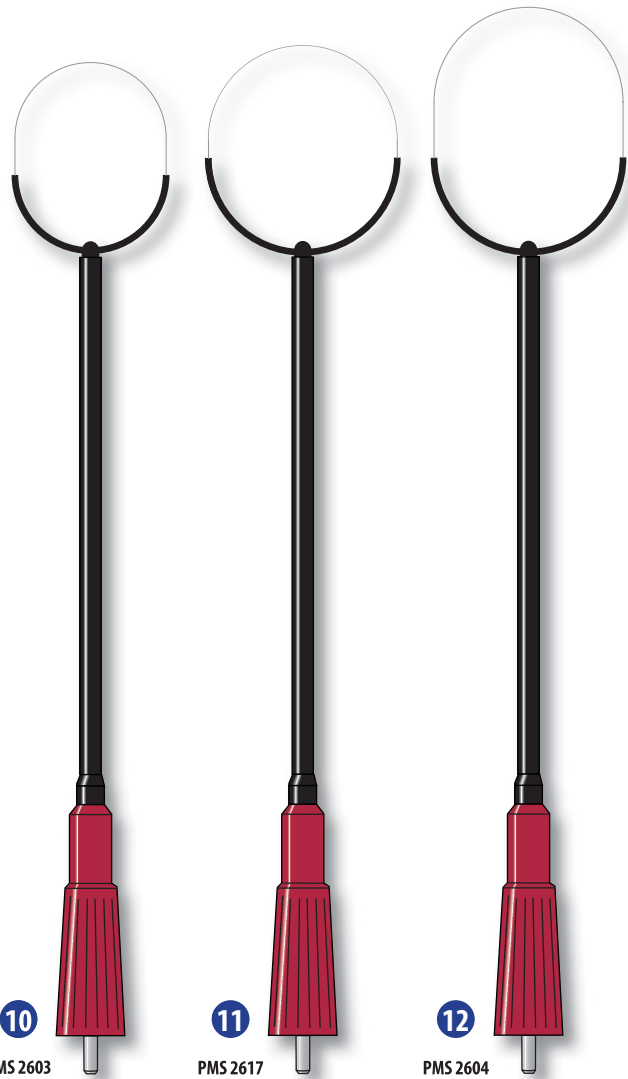

Loop shape may be adjusted

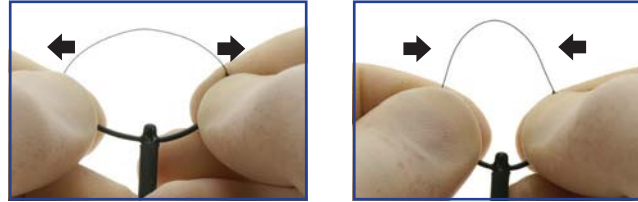

The Prima Flexiline Pro Lletz system has been designed for increased precision in the treatment of cervical dysplasia. This is achieved by improving accuracy of removal of biopsy or lesion from the transformation zone.

Reinforced coated metal frame permits shaping of the active tip allowing the surgeon to manipulate the loop to the specific requirements of the patient, minimising removal of good tissue.

The symmetrical shape of the loop allows for clear vision of the active tip when in procedure, resulting in a more accurate tissue dissection, especially when working with larger loops where over removal can cause problems with misshaping of the cervix.

Non-stick ball electrodes minimise adhesion and prevent bleeders from resuming when concluding procedure. This further allows for fulguration to be achieved with greater success rates than with standard stainless steel ball electrodes.

Safety shield



Safety shield designed to prevent potential arcing between the electrode and fingerswitch.

Sterile packed



Sterile packaged in protective electrode tray.


Single use non-stick ball electrodes



10	Loop electrode, 20 mm x 15 mm, 11 cm,	pack 10	PMS2603
11	Loop electrode, 25 mm x 15 mm, 11 cm,	pack 10	PMS2617
12	Loop electrode, 25 mm x 20 mm, 11 cm,	pack 10	PMS2604

FLEXILINE PRO SINGLE USE NON-STICK BALL ELECTRODES

13	3 mm ball electrode, 11 cm,	pack 10	PMS2628
14	4 mm ball electrode, 11 cm,	pack 10	PMS2627
15	5 mm ball electrode, 11 cm,	pack 10	PMS2626

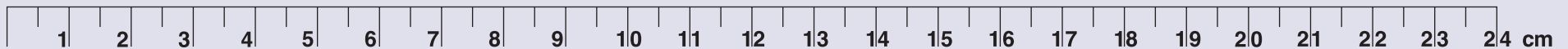

view our online catalogue at

www.prima-medical.com

Call today
to receive your
FREE
Prima catalogue

